


## APPLICATION NEWS

### Catecholamines and Related Compounds

Column: AquaSep 3 micron 100 x 4.6mmID  
Catalog #: 125121-AQS  
Mobile Phase: 50mM Phosphoric Acid/Methanol (93:7)  
Flow: 1.3mL/min.  
Detection: uv @ 225nm  
Injection Volume: 5uL

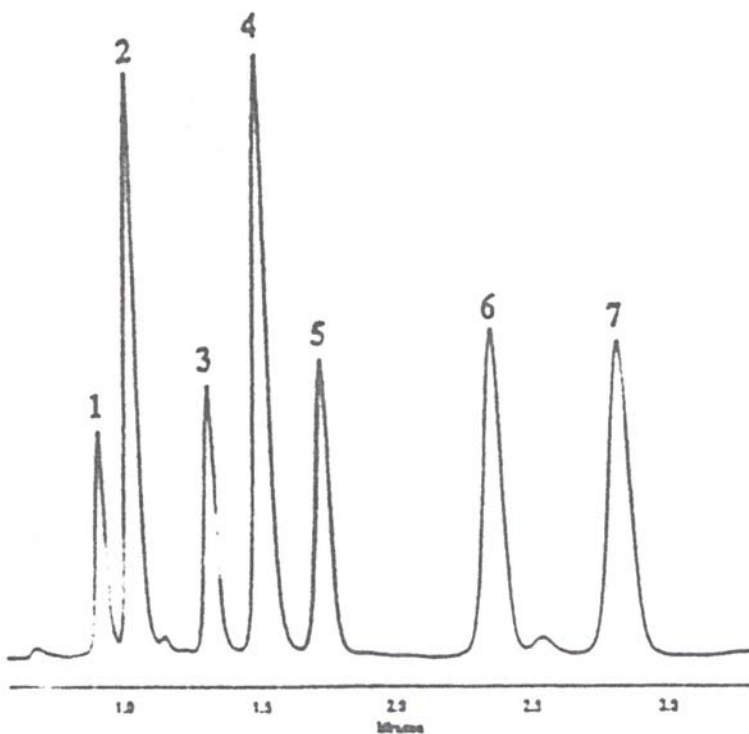

Sample: 1. Norepinephrine (8 ng)  
2. Epinephrine (20ng)  
3. 3-Hydroxylamine (8 ng)  
4. Metanephrine (19 ng)  
5. L-DOPA (9 ng)  
6. Tyrosine (9 ng)  
7. Serotonin (9 ng)